

Simply Slim Classic®



ReadyWear
Amplification
ITC

DEVICE FEATURES

- Small, in-the-ear design
- #10 battery
- Soft silicone body
- Class D Circuitry
- Hypo allergenic
- Volume Control

CONSUMER APPEAL

- Cosmetic Appeal
- Economical Price
- No Ear Impression
- Fits men and women

Soft, petite size with precision volume control.


Simply Slim Classic is specifically designed for a petite, or narrow ear canal shape. Simply Slim utilizes SoftEar Technology, a design that embodies sophisticated micro-electronic circuitry in a solid body of medical grade silicone. Due to the soft, flexible nature of the body, the device seats better in the ear and stays in place to minimize feedback and distortion. The silicone body is noticeably narrow in design when compared to the traditional Simply Soft device. It is 100% hypo-allergenic and constantly flexes with jaw movement. Simply Slim Classic offers a manual volume control allowing users to set an exact volume level.



DATA SPECIFICATIONS

Peak SSPL 90	111 dB SPL
HFA SSPL 90	110 dB SPL
Peak Full-on Gain	23 dB
HFA Full-on Gain	20 dB
Frequency Response.....	250-6900 Hz.
Total Harmonic Distortion.....	<3.5%
Battery Drain	0.21 mA
Battery Life (10A)	350 Hrs.
Equivalent Input Noise	<34 dB

TYPICAL FITTING RANGE

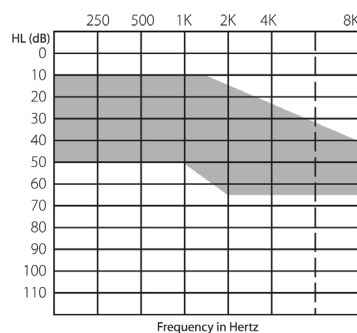

Simply Slim Classic®

Simply Slim Classic Data

The chart below reveal Simply Slim Classic smooth, wide-band frequency response, modest gain, and robust output. The smooth frequency response ensures clear, clean, high fidelity amplification. The robust output and modest gain combine to provide a pleasant sense of hearing enhancement without overpowering the hearing system or resulting in distortion due to sound saturation of the receiver.

Effects of volume control rotation demonstrates smooth, predictable gain adjustments provided by the Simply Slims manual volume control, allowing for benefit over a wide range of hearing losses.

