

Simply Soft Classic[®]

ReadyWear
Amplification

ITC

DEVICE FEATURES

- Small, in-the-ear design
- #10 battery
- Soft silicone body
- Class D Circuitry
- Hypo allergenic
- Volume Control

CONSUMER APPEAL

- Cosmetic Appeal
- Economical Price
- No Ear Impression
- Fits men and women



The original soft hearing device.


Simply Soft Classic, originally patented in the 1990's, is still the ultimate soft hearing device. Simply Soft was the first device to Implement SoftEar Technology, a design that embodies sophisticated micro- electronic circuitry in a solid body of medical grade silicone. Due to the soft, flexible nature of the body, the device seats better in the ear and stays in place to minimize feedback and distortion. The silicone is hypo allergenic and constantly flexes with jaw movement. Simply Soft Classic is equipped with high fidelity, high efficiency, Class D Circuitry.



DATA SPECIFICATIONS

Peak SSPL 90	111 dBSPL
HFA SSPL 90	110 dBSPL
Peak Full-on Gain	23 dB
HFA Full-on Gain	20 dB
Frequency Response.....	250-6900 Hz.
Total Harmonic Distortion.....	<3.5%
Battery Drain	0.21 mA
Battery Life (10A)	350 Hrs.
Equivalent Input Noise	<34 dB

TYPICAL FITTING RANGE

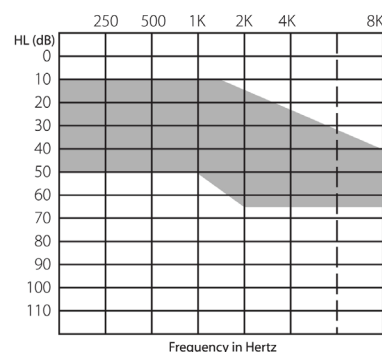

Where there's so much to hear about...

www.generalhearing.com

ghi General Hearing
Instruments, Inc.

P.O. Box 23748
New Orleans, LA 70183-0748
800.824-3021 (504) 733-3767

Simply Soft®

Simply Soft Classic Data

The chart below reveal Simply Soft Classic smooth, wide-band frequency response, modest gain, and robust output. The smooth frequency response ensures clear, clean, high fidelity amplification. The robust output and modest gain combine to provide a pleasant sense of hearing enhancement without overpowering the hearing system or resulting in distortion due to sound saturation of the receiver.

Effects of volume control rotation demonstrates smooth, predictable gain adjustments provided by the Simply Soft's manual volume control, allowing for benefit over a wide range of hearing losses.

