

Simply Tranquil™ I



The only all-soft, mini ITC tinnitus noise generator for extended wear patients.

The Simply Tranquil I is a ready-to-wear tinnitus noise generator device offering digital signal quality. Patients will appreciate same day fitting convenience - no ear impressions are needed!

Device Features

- Open-ear venting
- Rotary VC
- Low battery alert
- 0 dB SL floor
- #10 battery
- 200 hrs battery life
- No programming required
- Hypoallergenic soft body



Soothing sounds of the Tranquil™ brand noise generator that hearing professionals and patients have trusted since 1996.

The Simply Tranquil™ I is the preferred choice of tinnitus patients who need a noise generator they can wear all day and throughout the night! It is one of the smallest, most comfortable low-level sound generator on the market today. It offers natural sound quality due to the open canal format in a small ITC (In-The-Canal) design. The all-soft body will never harden over time and provides a comfortable, unoccluded fit.

The patented Tranquil circuitry provides the user with stable, broad-band noise controlled by a rotary volume control for accurate adjustments. Due to the soft, flexible nature of the body, the device seats comfortably in the ear and constantly flexes with jaw movement.

The Simply Tranquil I is powered by #10 battery offering approximately 200 hours of battery use. Be sure to specify preference of either general white noise or general pink noise when ordering the Simply Tranquil I.

Device Data

Both noise generator circuits provide the user with a 0 dB SL noise floor setting.

* Specify preference of General white or Pink noise sound generator at time of Tranquil device order.


GENERAL WHITE NOISE RESPONSE

Source.....OFF
 Peak.....72.9 dB SPL
 Peak Frequency.....300 Hz
 RMS OUT.....80.4 dB SPL
 Battery Drain.....0.40 mA
 Battery Life (10A).....200 Hrs.

GENERAL PINK NOISE RESPONSE

Source.....0.0 dB
 Peak.....79.7 dB SPL
 Peak Frequency.....300 Hz
 RMS OUT.....82.6 dB SPL
 Battery Drain.....0.40 mA
 Battery Life (10A).....200 Hrs.

FREQUENCY RESPONSE

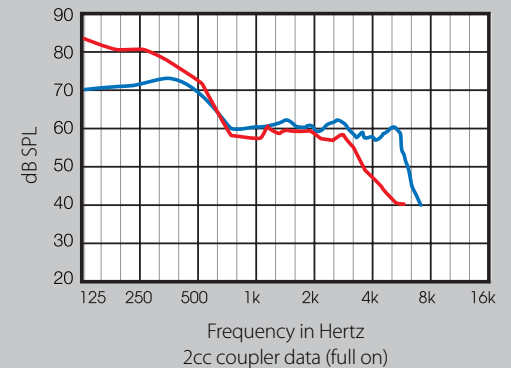

BLUE CURVE = General white noise

RED CURVE = General pink noise

ghi General Hearing
 Instruments, Inc.

175 Brookhollow Esplanade
 Harahan, LA 70123
 800.824.3021 504.733.3767

Hearing Made Simple
www.generalhearing.com

