

Tranquil™ I OTE

tinnitus noise generator



A mini OTE tinnitus sound generator that offers broad-band general or pink noise.

The Tranquil I OTE utilizes advanced digital signal processing with a choice of general white or general pink noise integrated into a cosmetically appealing mini OTE case design.

Device Features

- No programming required
- #10 battery
- Rotary volume control
- 200 hours battery life
- Low battery alert
- Open-ear dome
- 0 dBSL floor



Quality Tranquil™ brand noise generator in a mini case design with micro polytube and popular rotary volume control.

A robust volume wheel controls the Tranquil I which provides easy and accurate adjustment. The device offers a zero noise floor, which can be a critical feature in setting the volume of the device.

The Tranquil OTE (Over-The-Ear) device employs an ultra-cosmetic design combined with a micro polytube free-field sound dome which will enable same-day fitting scenarios. The open-ear design of the soft sound dome tip leaves the ear canal open, delivering clean stable sound while allowing the wearer to fully appreciate surrounding sounds in a normal manner. Be sure to specify preference of either general white or general pink noise when ordering the Tranquil I OTE.

The Tranquil I uses a #10 battery providing a typical 200 hours of use.

Device Options...

- Choice of case colors
- Tube length: 0-3 sizes available
- Dome size: small is standard, medium or large available
- Polytube fitting gauge at www.generalhearing.com



Device Data

Both therapy sounds will provide the user with a 0 dBSL noise floor setting.


*** Specify preference of General White or Pink noise sound generator at time of Tranquil device order.**

GENERAL WHITE NOISE RESPONSE

Source.....OFF
 Peak.....77.2 dB SPL
 Peak Frequency.....400 Hz
 RMS OUT.....88.3 dB SPL
 Battery Drain.....0.40 mA
 Battery Life (10A).....200 Hrs.

GENERAL PINK NOISE RESPONSE

Source.....OFF
 Peak.....82.7dB SPL
 Peak Frequency.....200 Hz
 RMS OUT.....86.9 dB SPL
 Battery Drain.....0.4 mA
 Battery Life (10A).....200 Hrs.

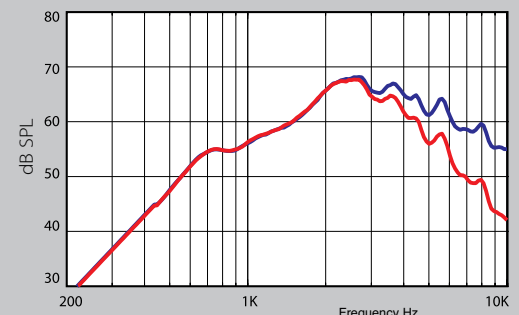

KEMAR Real-Ear (full on)

BLUE CURVE = General white noise

RED CURVE = General pink noise

ghi General Hearing
Instruments, Inc.

175 Brookhollow Esplanade
 Harahan, LA 70123
 800.824-3021 504.733.3767

Hearing Made Simple
www.generalhearing.com

