

# Tranquil™ Custom II



The unique Canal Open Ear (COE) tinnitus now with a multi-function switch.

The Tranquil II Custom's patented design fits comfortably, offers unoccluded natural sound quality, and allows users to switch between general white or pink noise.

### Device Features

- Open-ear venting
- Memory/volume selector
- Tranquil circuitry
- 0 dB SL noise floor
- Low battery alert
- #10 battery
- 200 hours battery life
- High output available
- Digital circuitry



### The updated version of the popular Tranquil COE design.

Just like the Tranquil™ I Custom, the Tranquil II is a canal-open-ear, tinnitus noise generator. The device offers advanced digital circuitry with the same quality broadband noise as featured throughout the Tranquil product line. Tranquil II is a custom-fit device that occupies only fifty percent of the ear canal. This unique shell design allows for a comfortable, unoccluded fit allowing patients who do not require amplification to hear in a normal manner.

Tranquil II offers a memory/volume toggle that is simple to use and combines two functions into one control. Holding the toggle in either direction will result in a smooth memory change; momentary "clicking" the toggle in an upward or downward motion results in discreet 1dB changes in volume.

The device provides a 0 dB SL noise floor volume control - a setting at which no audible noise is produced when volume is set at lowest setting. Tranquil II is powered by a zinc-air #10 battery that provides approximately 200 hours battery life.

### Device Data

Both therapy sounds provide the user with a 0 dB SL noise floor setting.


#### GENERAL WHITE NOISE RESPONSE

Source.....OFF  
 Peak.....72.9 dB SPL  
 Peak Frequency.....300 Hz  
 RMS OUT.....80.4 dB SPL  
 Battery Drain.....0.40 mA  
 Battery Life (10A).....200 Hrs.

#### GENERAL PINK NOISE RESPONSE

Source.....0.0 dB  
 Peak.....79.7 dB SPL  
 Peak Frequency.....300 Hz  
 RMS OUT.....82.6 dB SPL  
 Battery Drain.....0.40 mA  
 Battery Life (10A).....200 Hrs.

#### FREQUENCY RESPONSE

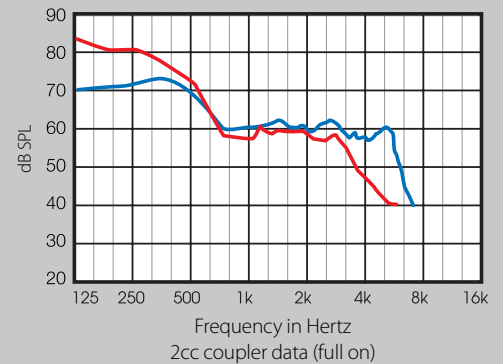

**BLUE CURVE** = General white noise  
**RED CURVE** = General pink noise



175 Brookhollow Esplanade  
 Harahan, LA 70123  
 800.824.3021 504.733.3767

Hearing Made Simple  
[www.generalhearing.com](http://www.generalhearing.com)

